


Fitting Instructions for 3M™ Filtering Facepiece Respirators

Issue Date 03/26/04


Wearing your filtering facepiece respirator

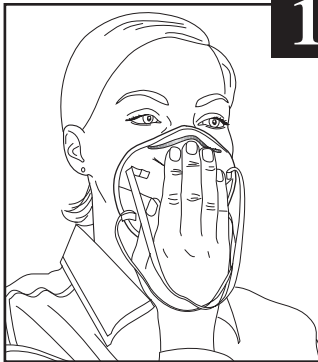


1 Place the respirator over your nose and mouth. Be sure the metal nose clip is on top. With model 8210, pre-stretch the straps before wearing.




2 Pull the top strap over your head until it rests on the crown of your head above your ears.

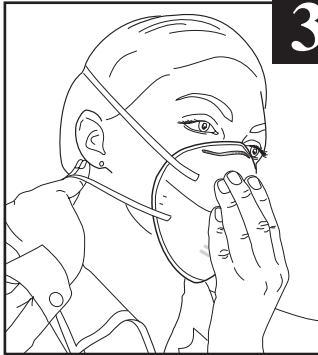

3 Pull the bottom strap over your head until it rests just below your ears.

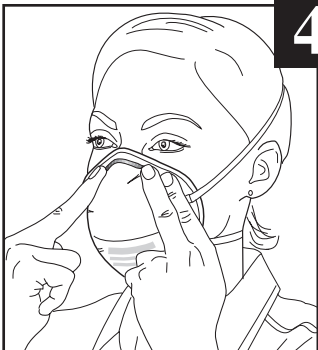

4 Using both hands starting at the top, mold the metal nose clip around your nose to achieve a secure seal. **Note:** Pinching the nosepiece using one hand may result in improper fit and less effective respirator performance. Use two hands.

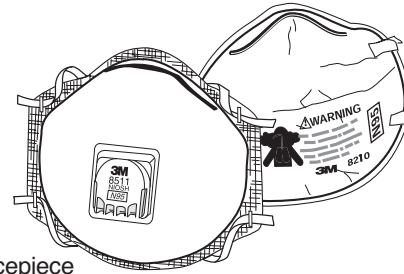


Filtering facepiece non-valved respirator


Filtering facepiece valved respirator

Check the seal of your filtering facepiece respirator each time you don the respirator.

1a For non-valved respirators Place both hands completely over the respirator and *exhale*. The respirator should bulge slightly. If air leaks between the face and face seal of the respirator, reposition it and readjust the nose clip for a more secure seal. If you cannot achieve a proper seal, **do not** enter the contaminated area. See your supervisor.

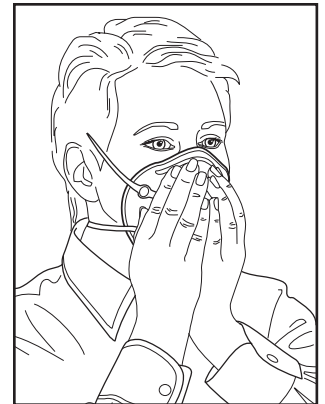

1b For valved respirators Place both hands over the respirator and *inhale* sharply. The respirator should collapse slightly. If air leaks between the face and face seal of the respirator, reposition it and readjust the nose clip for a more secure seal. If you cannot achieve a proper seal, **do not** enter the contaminated area. See your supervisor.

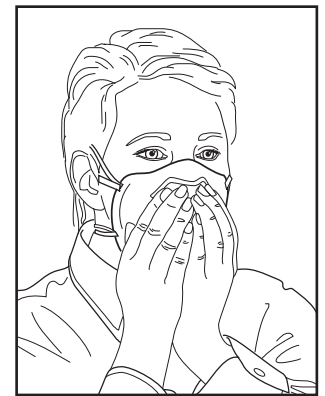



⚠ WARNING

These respirators help reduce exposure to certain particles. **Misuse may result in sickness or death.** Before use, the wearer must read and understand **User Instructions** provided as a part of product packaging. Time use limitations may apply. For proper use, see package instructions, supervisor or call 3M OH&ESD Technical Service in U.S.A., 1-800-243-4630. In Canada, call 1-800-267-4414.

Important

Before using these respirators, you must determine the following:

1. The type of contaminant(s) for which the respirator is being selected.
2. The concentration level of contaminant(s).
3. Whether the respirator can be properly fitted on the wearer's face. Do not use with beards, on other facial hair, or other conditions that prevent a good seal between the face and the faceseal of the respirator.
4. Before use of these respirators, a written respiratory protection program must be implemented, meeting all the requirements of OSHA 29 CFR 1910.134, including training, medical evaluation and fit testing.

For more information, please contact:

3M Occupational Health and Environmental Safety Division (OH&ESD)

In the U.S., contact:

Sales Assistance

1-800-896-4223

Technical Assistance

1-800-243-4630

Fax On Demand

1-800-646-1655

Internet

<http://www.3M.com/occsafety>

For other 3M products

1-800-3M HELPS

In Canada, contact:

3M Canada Company, OH&ESD

P.O. Box 5757

London, Ontario N6A 4T1

Sales Assistance

1-800-265-1840, ext. 6137

Technical Assistance (Canada only)

1-800-267-4414

Fax On Demand

1-800-646-1655

Internet

<http://www.3M.com/CA/occsafety>

Technical Assistance In Mexico

01-800-712-0646

5270-2255, 5270-2119 (Mexico City only)

Technical Assistance In Brazil

0800-132333

Fax On Demand O.U.S. Locations

1-651-732-6530



General Offices

3M Center

St. Paul, MN 55144-1000